


Shucman ²

Fibre Optic Laryngoscope Blades

*compatible with
conventional laryngoscope handles*

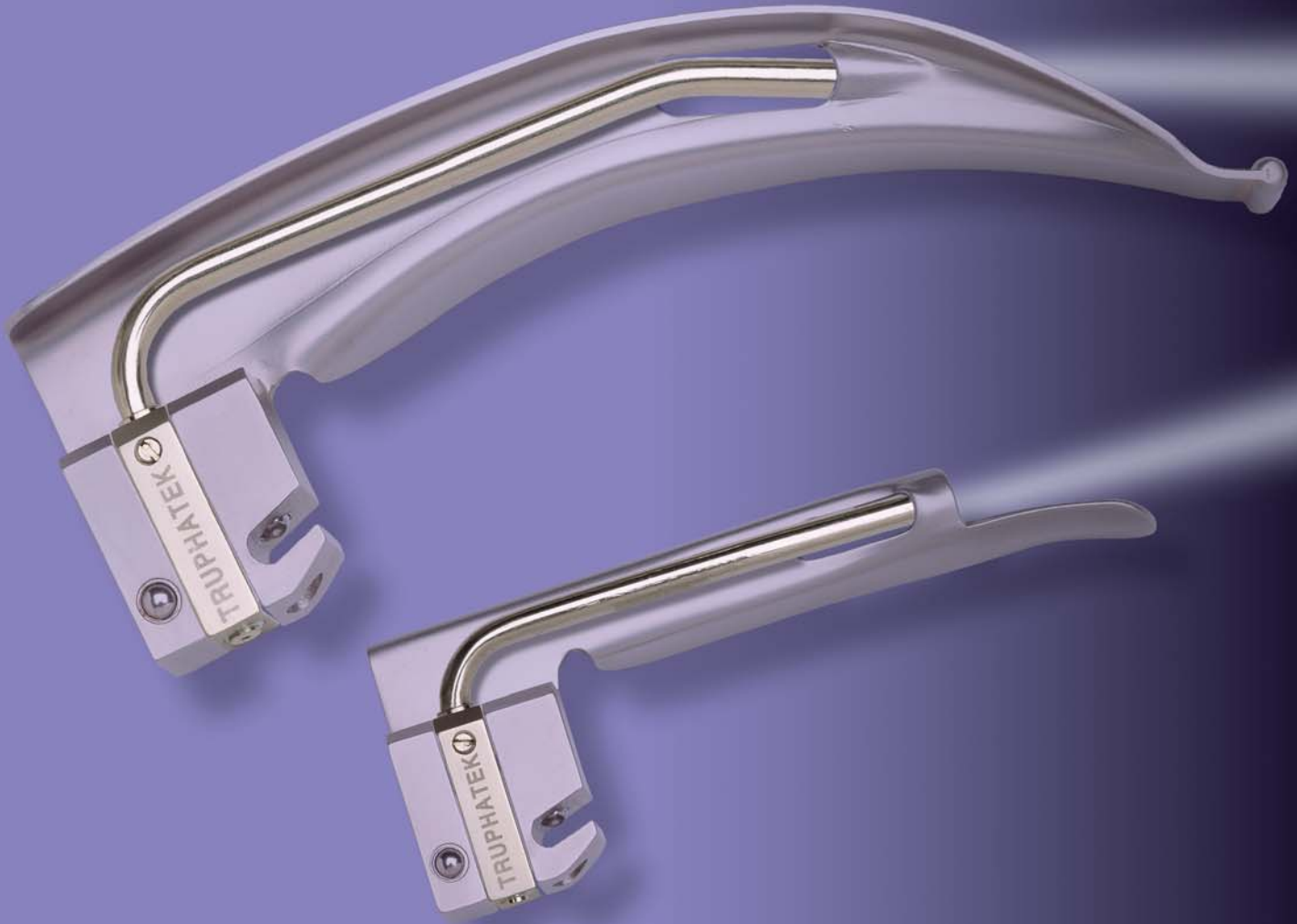

*Your easy “upgrade”
to Fibre Optics*

Conforming to ISO 7376-1 & ASTM F 965-99
Patents registered and pending worldwide
Approvals - ISO 9001, ISO 13485
FDA registered manufacturer



 **TRUPHATEK[®]**

5 Reasons to Use

3 Times More Light Output

1



Shucman™

- Light intensity in excess of 2,400 Lux
- Cold, white light illumination
- Halogen bulb located in base of the blade
- No patient risks

Warm light Blade

- Light intensity less than 750 Lux
- Warm pink bulb illumination
- Incandescent bulb at tip of blade
- Hot or broken bulbs put your patient at risk

Quality + Value + Customer service

2



Shucman™

- Leading world supplier using highest grade materials
- Each blade bears batch number
- Designed to be fully autoclavable
- Local service centre available
- Each blade carries a 5 year manufacturers' guarantee

Conforms to ISO 7376(1)

3



Shucman™

- All products are FDA and CE listed
- Each stainless steel blade conforms strictly to all international standards
- Can be used with all handles that conform to ISO 7376(1)

Upgrade Today

Unique Patented Replaceable FibreClip™



Shucman™

- Each FibreClip™ contains halogen bulb and fibre optic light guide and has a working life of at least 15 months
- If required FibreClip™ may be removed for cleaning, autoclaving and replacement
- FibreClip™ extends the blade's life indefinitely

Shucman™ Blade and **Dolphin™** Handle
Set variations available



Unique Dolphin™ Waterproof Handle



Dolphin™

- Allows immersion in cleaning solutions without removal of batteries
- With removal of batteries, Dolphin™ handles are fully autoclavable
- Long life with the highest quality

Technical Specifications



Size	Use	Overall Length (mm)	Distal Width (mm)	Order No.
Mac0	Neo-Natal	77	9	5240
Mac1	Neo-Natal	90	9	5241
Mac2	Infants	100	12	5242
Mac3	Adults	130	15	5243
Mac4	Adults	155	15	5244
Mac5	Adults	172	15	5255



Size	Use	Overall Length (mm)	Distal Width (mm)	Order No.
Miller00	Neo-Natal	64	11.5	5247
Miller0	Neo-Natal	75	11.5	5250
Miller1	Infants	102	11.5	5251
Miller2	Child	155	14	5252
Miller3	Adults	195	14	5253
Miller4	Adults	205	14	5254



Size	Use	Overall Length (mm)	Distal Width (mm)	Order No.
OxyMiller00	Neo-Natal	76	11.5	5248
OxyMiller0	Neo-Natal	87	11.5	5249



Sets Include *		
Handle	Blades with FiberClip™ Light Guide	Order No.
Dolphin™ Medium	Mac 1,2 & 3	5440 D
Dolphin™ Medium	Mac 1,2,3 & 4	5441 D
Dolphin™ Medium	Mac 2,3 & 4	5439 D
Dolphin™ Mini	Miller 00,0 & 1	5430 D
Dolphin™ Medium	Miller 2,3 & 4	5456 D

* Additional combinations available upon request.

Distributed by



Size	Use	Overall Length (mm)	Order No. Pack of 2 FibreClip™ Light Guide
Mac0	Neo-Natal	54	524002
Mac1	Neo-Natal	54	524102
Mac2	Infants	66	524202
Mac3	Adults	90	524302
Mac4	Adults	116	524402
Mac5	Adults	116	525502
Miller00	Neo-Natal	36	524702
Miller0	Neo-Natal	42	525002
Miller1	Infants	70	525102
Miller2	Child	109	525202
Miller3	Adults	146	525302
Miller4	Adults	156	525402

Size	Finish	Batteries	Code
UltraSafe™ Stubby	Knurled	AA	5523
Dolphin™ Medium	Knurled	C	5571
Dolphin™ Mini	Knurled	AA	5572

From Truphatek's® Warranty Statement:

All Truphatek® products are unconditionally guaranteed to be of good quality and free from manufacturing defects in material and workmanship for a period of 5 years from date of manufacture.

FibreClip™ fibre optic light guides are guaranteed for 15 months use from date of manufacture, as it appears embossed on the base of the light guide.

Replacement product will be provided without charge for any Truphatek® product not meeting this warranty."

Full text of Truphatek's® warranty can be viewed at our website www.truphatek.com or by request from Customer Services at helpdsk@truphatek.com

Shucman 2™, UltraSafe™, Dolphin™ and FibreClip™ are trademarks of Truphatek International and should not be used without prior written permission from Truphatek International Ltd., Israel.



Truphatek International Ltd. P.O.Box 8051 Netanya 42504, Israel
Tel: 972-9-8851155, Fax: 972-9-8851212, helpdsk@truphatek.com

© 2003-2007 TRUPHATEK® All rights reserved

SHUC2-2-05/07