



Long Life Power Handle

With Dependable, Consistent LED Light Output



INDESTRUCTIBLY ENGINEERED


- Stainless steel 316 prevents normal wear & tear caused by repeated mounting of blades
- Resistant to rust and chemicals
- Can be sterilized & autoclaved repeatedly by all prevalent systems

UNIQUE ERGONOMIC DESIGN

- One ideal size suitable for all patient anatomies
- Silicon cushioned grip

TESTED AND APPROVED FOR HIGHEST "UL" SAFETY REQUIREMENTS

- Unique controller protects battery life and charge efficiency
- Multi-coloured LEDs indicate battery integrity and charge status
- Classified for non-interference with devices in proximity

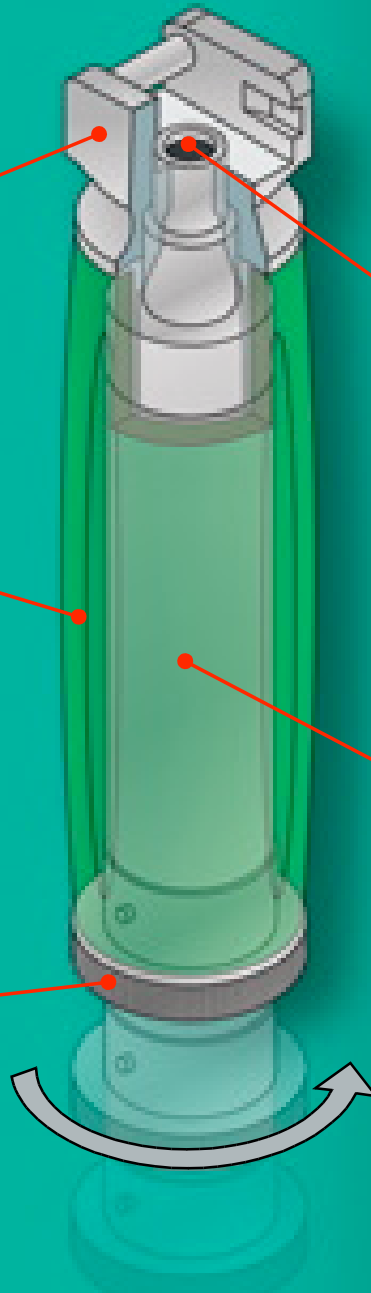

RELIABLE LIGHT SOURCE

- LED provides optimal bright bluish white light spectrum
- Unique, patented configuration emits uniform and focused light.
- Light output substantially unvaried over 5 year life span

RECHARGEABLE LITHIUM ION BATTERY

- No capacity deterioration or memory build-up
- Power for over 10 hours of continuous use between charges
- Recharge time equals use time
- RoHS compliant for banned materials

Compatible with all Green System laryngoscope blades
 Conforms to ISO 7376-3 and ASTM F 1195-99
 Approvals - ISO 9001, ISO 13485
 FDA registered manufacturer



Frequently Asked Questions

- **How many hours can I use the Tru-LED™ without recharging the battery?**
No less than 10 hours of continuous use.
- **What is the life span of the LED?**
Under optimal conditions, an LED has over 50,000 hours of use. Truphatek guarantees the LED for 5 years from date of manufacture.
- **What is the life span of the Lithium Ion Battery?**
Over 300 full cycles. One cycle is from a completely empty battery to a fully charged battery.
In addition, the battery is guaranteed for 5 years (see Truphatek's Warranty Statement at www.truphatek.com).
- **By using LED as the light source instead of Xenon bulb, what has my hospital saved?**
The cost of dozens of expensive Xenon bulbs. And, in addition, you'll be saving all the fuss and bother of replacing bulbs.
- **By using rechargeable batteries instead of alkaline batteries, what has my hospital saved?**
An average hospital disposes of hundreds of batteries which you will now be saving. And besides, you'll be protecting our environment by decreasing the number of batteries to be disposed.
- **How can I clean the Tru-LED handle?**
The outer sleeve of the handle is made of either aluminum or silicon padding. When the outer sleeve is separated from the inner cartridge, it is completely autoclave-able.
The inner cartridge containing the rechargeable battery should be wiped down gently with a dry cloth, taking care not to scratch the glass protecting the LED at the top of the cartridge.
- **Is the Tru-LED handle available in additional sizes?**
Yes, a distinctive, rechargeable handle is also available in a specially designed Paediatric size, intended primarily for use with infants. More details on www.truphatek.com.

Technical Specifications

Ordering Code	Product	Description
4220 E,A,B or U/J(*)	Tru-LED™ Handle set	Tru-LED™ ergonomical handle with rechargeable battery cartridge and power supply in clear plastic case
42143000	Spare Battery Cartridge	For Standard size handle
421460 E,A,B or U/J (*)	Spare Power Supply	Suitable for standard and infant size handles
Handle and Blade sets Upon request	Sets containing Macintosh or Miller Green Spec™ blades, with one Tru-LED™ handle, battery cartridge and power supply	Sets with 3 or 4 blades Sizes upon request



(*)E = European plug, A = Australian plug
B = British plug, U/J = USA and Japanese plug

From Truphatek's® Warranty Statement:

"All Truphatek® products are unconditionally guaranteed to be of good quality and free from manufacturing defects in material and workmanship for a period of 5 years from date of manufacture. Replacement product will be provided without charge for any Truphatek® product not meeting this warranty." Full text of Truphatek's® warranty can be viewed at our website www.truphatek.com or by request from Customer Services at helpdesk@truphatek.com
Tru-LED™, Truphatek® and GreenSpec™ are trademarks of Truphatek International and should not be used without prior written permission from Truphatek International Ltd., Israel.

Distributed by

 **TRUPHATEK®**
www.truphatek.com

Truphatek International Ltd. P.O.Box 8051 Netanya 42504, Israel
Tel: +972-9-8851155, Fax: +972-9-8851212, helpdesk@truphatek.com

© 2008 TRUPHATEK®

All rights reserved

TL-01-01/08